# armfield

# Hydraulics & Hydrology - S series



# Standard teaching and research flume – S6-MKIII

armfield

The Armfield S6-MKIII laboratory flow channel is one of the most important tools available to the hydraulics or civil engineer whether engaged in teaching basic principles or researching solutions to practical problems.

Many applications in fluid mechanics are associated with the flow of water through an open channel where the water has a free surface that is exposed to the air at atmospheric pressure.

The flumes are available in different lengths from 5 to 17.5 meters increasing in 2.5 meter increments. Armfield flumes are installed in educational and research establishments throughout the world.

TILTING UP TO 17.5 METERS MODULAR DESIGN CONTROL AND DATA LOGGING VIA PC SEDIMENT TRANSPORT OPTIONS MANUAL OR ELECTRICAL JACKING





#### Features / benefits

- Accurate for education and research
- Designed for ease of visibility:- toughened glass sides, slimline side supports and comfortable viewing height

UK office - email: sales@armfield.co.uk tel: +44 (0) 1425 478781 (for ROW) USA office - email: info@armfield.inc tel: +1 (609) 208-2800 (USA only)

- Modular construction supplied in pre-glazed sections for rapid and easy site assembly
- Wave generation options which can be used to propagate random or regular waves in the working section

Issue: 2	Applic	ations
URL: http://www.armfield.co.uk/flumes	ME	CE
/e reserve the right to amend these specifications without prior notice. E&0E © 2021 Armfield Lt	d. All Rights	s Reserved

# armfield.co.uk


### S6-MK-III Glass Sided Tilting Flumes

#### Description

There are numerous design features associated with Armfield flumes \* Not all features are appropriate for every channel.

- Accurate for education and research
- Extremely stable design, no user adjustments required to the flume bed
- Floor space requirements reduced to a minimum
- Fabricated high precision stainless steel channel bed
- Quick conversion to closed-loop recirculation for sediment transport studies
- Precision screw jacks provide accurate slope adjustment with minimum effort (powered jacks available as an option)
- Adjustable instrument rails with positioning scales fitted over the whole working length
- Fully profiled inlet tank fitted with stilling and smoothing devices

- Discharge tank with adjustable overshot weir and draft tube to avoid splashing and enhance noise reduction
- Modular construction supplied in pre-glazed sections for rapid and easy assembly on site
- Wave generation options, both regular and random
- Comprehensive range of optional accessories, instruments and models available
- Non-corroding durable GRP tanks throughout
- Transparent sides are of toughened glass, which is extremely strong, abrasion resistant, dimensionally stable, does not discolour or scratch and is inherently safe
- Working section allows adjust-ability, enabling extremely accurate setting
- Under frame designed to reduce load deflections to a minimum
- Close tolerances specified and achieved.

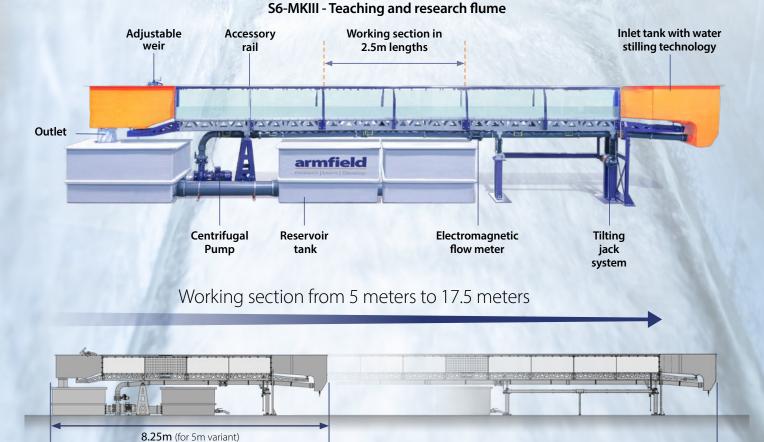

20.75m (for full 17.5m variant)



The most important aspect of a tilting flume is retaining the integrity of the working section.

To achieve this requires an extremely rigid design which ensures almost no deflection regardless of load or tilt.




#### **Experiments and Research**

The Armfield S6-MKII flume has been developed during 30 years of continuous production, and examples are installed in educational and research establishments throughout the world.

The flumes are available in different lengths to suit the application, short versions for basic investigations and longer versions for investigations of gradually varied flow profiles with non-uniform channel flow.



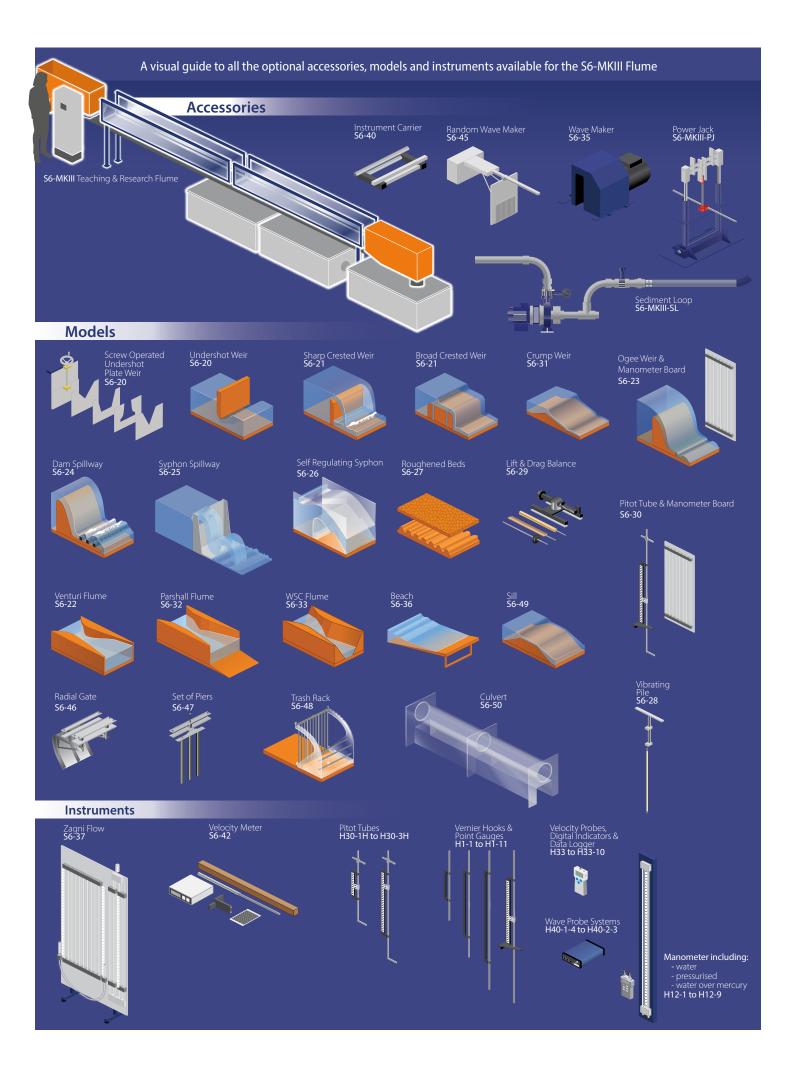

### **Ordering specification**

- ► A self-contained glass sided tilting flume for fluid mechanics laboratory experiments, project work and research activities
- The flume working channel is assembled from modular sections of 2.5m length. A wide choice of standard lengths are available from 5m upwards
- The flume cross-section is 300mm wide by 450mm deep ►
- A fabricated high precision stainless steel bed provides excellent strength and rigidity, eliminating the need for a separate underframe. No adjustments other than the jacking stations are necessary in order to set up and maintain the equipment, achieving typical bed deformations better than 1mm
- Each flume incorporates a discharge tank fitted with an adjustable overshot weir and draught tube to avoid splashing and noise
- An electro-magnetic flow meter is incorporated as standard
- A comprehensive range of optional accessories and instruments is ► available to supplement the capabilities of the basic flume
- Closed-loop recirculation is an option for sediment transport studies

# **Technical specifications**

Walls	Toughened glass		
Bed	Exclusively fabricated from stainless steel		
End tanks	GRP (Glass Reinforced Plastic)		
Sump tanks & pipework	PVC (Polyvinyl chloride) & PE (polyethylene)		
Pump	Close-coupled centrifugal		
Flow regulation valve	Hand wheel operated butterfly		
Flow meter	Electro-magnetic		
Maximum flow rate	38 Litres/sec		
Bed stability	1.0mm (typical) at 400mm water depth		
Side wall stability	0.8mm (typical) at 400mm water depth		
Width	0.3m		
Depth	0.45m		

#### Requirements

#### Scale

- لې کې SPh
- Electrical supply 3Ph, 50-60Hz
- ► The user must have access to a PC
- Software requires the user to have a PC running Windows 7 ► or above with a spare USB port

### **Tilting configurable modular flumes S6-MKIII**

#### Working section dimensions

Туре	Tilting or standard flume		
Width	0.3m		
Depth	0.45m		
Length (in 2.5m modular increments)	5m - 17.5m		

#### Packed and crated shipping specifications

Model	Volume	Gross weight	
S6-MKIII-5M	18m <sup>3</sup>	2600Kg	
S6-MKIII-7.5M	22m <sup>3</sup>	2900Kg	
S6-MKIII-10M	27m <sup>3</sup>	3200Kg	
S6-MKIII-12.5M	29m <sup>3</sup>	4400Kg	
S6-MKIII-15M	31m <sup>3</sup>	4700Kg	
S6-MKIII-17.5M	33m <sup>3</sup>	5000Kg	

# **Experimental Models & Instrumentation**

A comprehensive range of optional accessories, models and measuring instruments are available for selection. These provide the basis for a large number of practical experiments in open channel flow including the use and operation of regulating and gauging structures. Non-corroding materials have been used to reduce maintenance time and increase the working life of the models.

- ► S6-20: Plate Weirs
- ► S6-21: Broad Crested Weirs
- ► S6-22: Venturi Flume
- S6-23: Ogee Weir & Manometer Board
- ► S6-24: Dam Spillway Models
- ► S6-25: Syphon Spillway
- ► **S6-26:** Self-regulating Syphon
- S6-27: Roughened Beds
- S6-28: Vibrating Pile ► S6-29: Lift & Drag Balance & Models
- S6-30: Pitot Tube & Manometer Board
- ▶ S6-31: Crump Weir

- ► S6-32: Parshall Flume S6-33: WSC Flume
- S6-35: Wave Generator
- S6-36: Beach
- S6-37: Zagni Flow Monitoring Systems
- ► S6-40: Instrument Carrier S6-42: Velocity Meter
- and Mountings
- S6-45: Random Wave Maker
- S6-46: Radial Gate
- S6-47: Set of Piers
- S6-48: Trash Rack
- S6-49: Sill
- S6-50: Culvert
- ▶ S6-MKIII-DTA-ASUITE Software Control and Data Acquisition Package
- ► (S6-MkIII-DTA-ALITE) Data Logging and Instrumentation System

#### **Tilt parameters**

•						
Working Section	5M	7.5M	10M	12.5m	15m	17.5m
(+)%	5	5	4.5	3.6	2.9	2.5
(-)%	2.1	1.4	0.9	0.7	0.5	0.5
Total	7.1	6.4	5.4	4.3	3.4	3

### **Ordering codes**

#### S6 MKIII flumes - cross section 300mm wide x 450mm deep

S6-MKIII-5M-C S6-MKIII-5M-D S6-MKIII-7.5M-C S6-MKIII-7.5M-D S6-MKIII-10M-C S6-MKIII-10M-D S6-MKIII-12.5M-C S6-MKIII-12.5M-D S6-MKIII-15M-C S6-MKIII-15M-D S6-MKIII-17.5M-C S6-MKIII-17.5M-D S6-MKIII-SL S6-MKIII-PJ S6-MKIII-DTA-ASUITE S6-MKIII-DTA-ALITE

Self Contained 5mtr Flume 415V/3Ph/50Hz\* Self Contained 5 mtr Flume 208V/3Ph/60Hz\* Self Contained 7.5mtr Flume 415V/3Ph/50Hz\* Self Contained 7.5 mtr Flume 208V/3Ph/60Hz\* Self Contained 10 mtr Flume 415V/3Ph/50Hz\* Self Contained 10 mtr Flume 208V/3Ph/60Hz\* Self Contained 12.5mtr Flume 415V/3Ph/50Hz\* Self Contained 12.5 mtr Flume 208V/3Ph/60Hz\* Self Contained 15 mtr Flume 415V/3Ph/50Hz\* Self Contained 15 mtr Flume 208V/3Ph/60Hz\* Self Contained 17.5mtr Flume 415V/3Ph/50Hz\* Self Contained 17.5 mtr Flume 208V/3Ph/60Hz\* Sediment Loop for S6-MKIII flume all lengths Power Jacks for S6-MKIII flume all lengths Software Control and Data Acquisition Package Data Logging & Instrumentation Package \* includes Manual Jacks, Control Console with Pump and storage tanks

## S6-MKIII standard warranty applies with this product



# Aftercare

Installation Commissioning Training Service and maintenance Support: armfieldassist.com

# **Knowledge** base

> 28 years' expertise in research & development technology > 50 years' providing engaging engineering teaching equipment Benefit from our experience, just call or email to discuss your laboratory needs, latest project or application.